



Free-standing room covering rail system

– Easy, practical and safe

- **Easy to install and remove**
- **Suitable for temporary lifting needs**
- **Suitable for reuse**
- **Height-adjustable – max. height 284 cm (9'3")**
- **Requires no structural changes to the room or installation in walls or ceiling**
- **Frees up floor space and allows more room for the user and carer**

A mobile rail system that can be installed and removed quickly and easily as well as moved from place to place.

The system requires no structural changes to the room or installation in walls or ceiling, making it well-suited to temporary lifting needs in private homes, for example.

The free-standing rail system can easily be installed and removed by 1-2 people without the need for special tools.

The rail system is available in two standard dimensions, 3 x 3 meters (approx. 10' x 10') and 4 x 4 meters (approx. 13' x 13'), and is height adjustable up to a maximum of 284 cm (9'3").



Guldmann Inc.
 14401 McCormick Drive, Unit A
 Tampa, FL 33626
 Tel. 800 664 8834
 Tel. 813 880 0619
 Fax 813 880 9558
 info@guldmann.net
 www.guldmann.net

Free-standing room covering rail system, height adjustable

Technical specifications

Functions

Max. capacity : 770 lbs (350 kg)

Weight and Materials

Free-standing room covering

rail system, 3 x 3 m : 332 lbs (151 kg)
 4 Vertical poles : 24 lbs (11 kg), each
 4 B-rails : 36 lbs (16,4 kg), each
 4 top corner brackets : 11 lbs (4,8 kg), each
 1 Traverse rail
 with trolleys : 49 lbs (22,2 kg)

Free-standing room covering

rail system, 4 x 4 m : 392 lbs (178 kg)
 4 Vertical poles : 24 lbs (11 kg), each
 4 B-rails : 48 lbs (21,8 kg), each
 4 top corner brackets : 11 lbs (4,8 kg), each
 1 Traverse rail
 with trolleys : 61 lbs (27,6 kg)

Rails and poles : Anodized aluminium
 Brackets : Powder varnished steel

Measurements

Height min. : 2040 mm
 Height max. : 2840 mm

Length, 3 x 3 m : 3200 mm
 Width, 3 x 3 m : 3200 mm




Length, 4 x 4 m : 4200 mm
 Width, 4 x 4 m : 4200 mm

Labelling

The product is manufactured in compliance with the Council Directive 93/42/EEC of June 14th 1993, including amendments, as medical device class 1.

Item Number/Model

555262 : Free-standing room covering rail system, height adjustable – 3 x 3 m
 555570 : Free-standing room covering rail system, height adjustable – 4 x 4 m

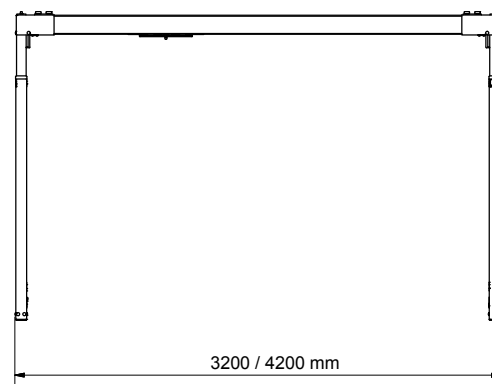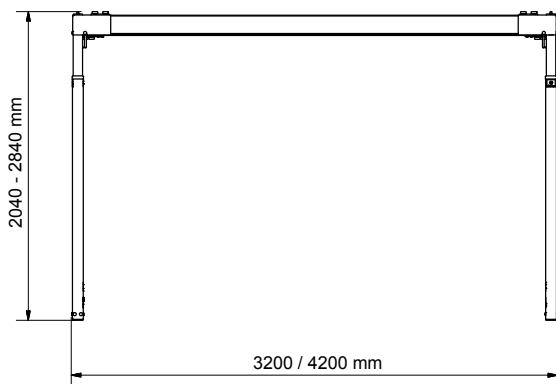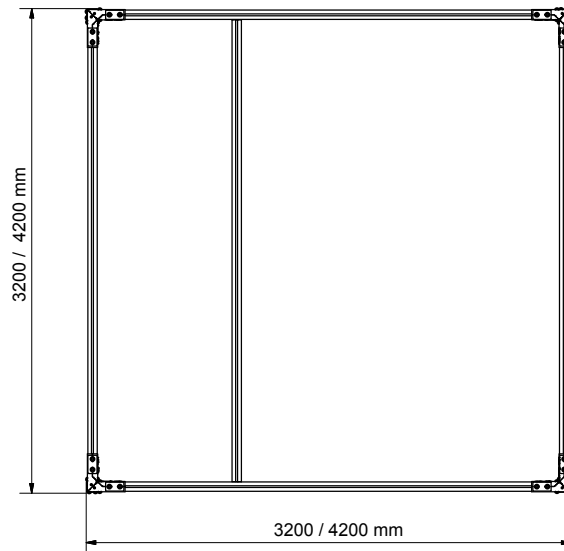

Time to care